




Brochure

Nexis[®] MIS Ø 2.7 mm: beveled compressive screws

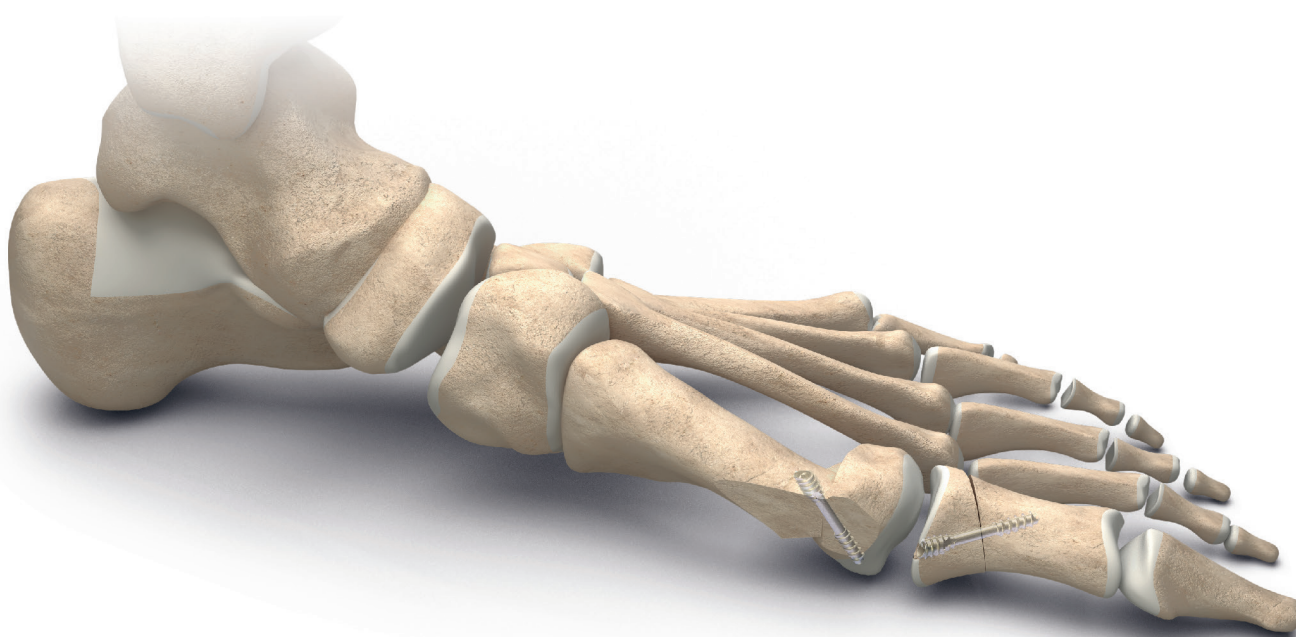


- . Specific & universal Exact2-T Recess
- . Elliptic beveled head
- . *Extra-Sharp* Ø 1.2 K-wire

The Nexis[®] MIS Ø 2.7 mm is a titanium alloy **cannulated compressive screw** featuring an **elliptic beveled head**.


Keys features

1 Cannulated screw: *Extra-Sharp Ø 1.2 K-Wire*


2 Elliptic beveled head

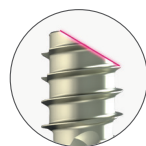


Bevel (30°):
Maximization of cortical anchorage and preservation of soft tissue

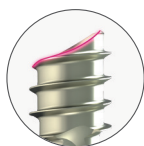


Elliptic bevel:
Allows for additional angular rotation to preserve the burial of the head

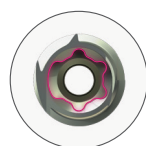

Exact2-T Recess:
Specific & universal

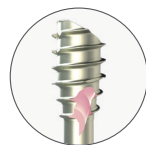


4 Positive locking channels

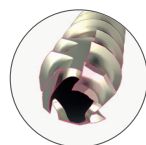

5 Self-drilling & self-tapping:
. Penetrating sharp tips
. Facilitated insertion

6 Deep dual thread

. Maximized anchorage & compression

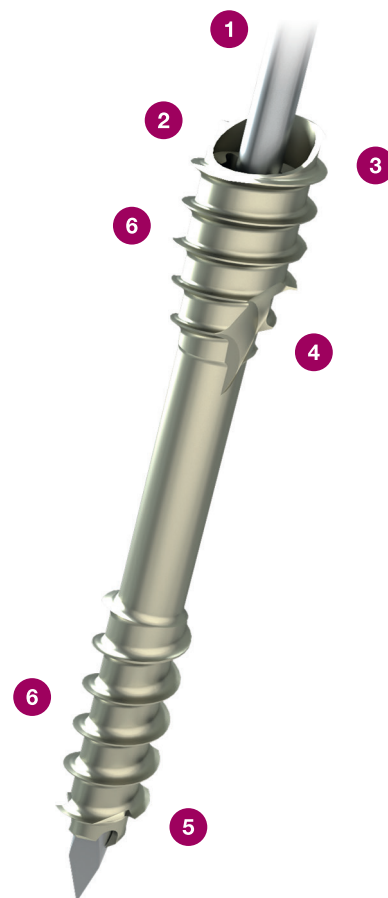

Optimized range
Screws available in lengths
of 14 to 30 mm

Optimized instrumentation

Extra-Sharp Ø 1.2 K-Wire

Ø 1.2 mm: provides greater stiffness*
Sharpened tip: allows for easy angulated insertion




Exact2-T8 screwdriver tip

Exact2-T8 recess:
. **Specific:** easy indexing of the Exact2-T8 AO screwdriver tip
. **Universal:** possible removal with standard instrumentation



Exact2-T8 recess ensures the unique positioning of the screwdriver tip.

The laser marking and bevel of the Exact2-T8 AO screwdriver tip facilitate the proper placement of the screw head.

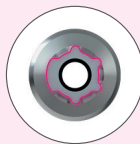


Laser marking

Bevel

*Greater stiffness than a Ø0.9 k-wire.



Nexis® MIS Mini Tray

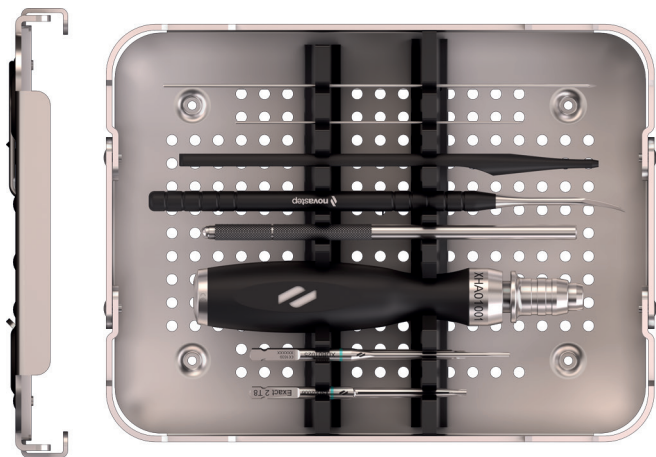

Extra-sharp K-wire & Exact2-T8 AO screwdriver tip: optional in the Nexis® ForefootCOMPLETE & ForefootEXACT trays

References

Nexis® MIS Ø 2.7 mm screws



Length (mm)	Reference
14 mm	SC090014
16 mm	SC090016
18 mm	SC090018
20 mm	SC090020
22 mm	SC090022
24 mm	SC090024
26 mm	SC090026
28 mm	SC090028
30 mm	SC090030


K-Wire Ø 1.2 mm

Description	Reference
K-Wire Ø 1.2 lg 80 TR/RD *	33-T10-R-12-080 *
K-Wire Ø 1.2 lg 100 TR/RD	33-T10-R-12-100
K-Wire Ø 1.2 lg 150 TR/RD **	33-T10-R-12-150 **
* Optional in the ForefootCOMPLETE	
** Optional in the Nexis® MIS Mini tray	


Indications

The osteosynthesis screws are indicated for arthrosis, hallux valgus and other bone alignment defaults (pes cavus, flatfoot, calcaneal valgus, malalignment secondary to previous trauma).

Chevron Osteotomy

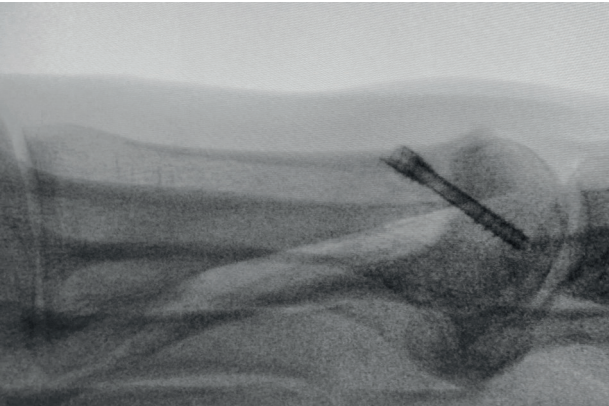

Akin Osteotomy

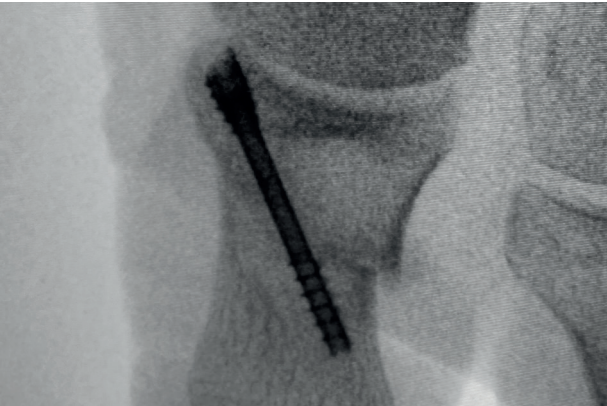

Please Note:
Carefully read the enclosed Instructions For Use (IFU) and all packaging label information. Devices: Implants: Class IIb-CE1639 / Instruments: Class I / Class IIa-CE1639.

Novastep:
2, Allée Jacques Frimot - 35000 RENNES - France
Tel: + 33 (0) 2 99 33 86 50 / Fax: + 33 (0) 9 70 29 18 95
contact@novastep-ortho.com / www.int.novastep.life

Reference: NEX-MIS-L-Ed2-10-23-EN