



FEATURES
BROCHURE

enovis™

PECA[®] COMPRESSIVE

Ø4 BEVELED COMPRESSIVE IMPLANTS

SCREW SYSTEMS



MULTIPLE INDICATIONS

The osteosynthesis screws are indicated for arthrosis, hallux valgus, metatarsalgia, and other bone alignment defaults (pes cavus, flatfoot, malalignment secondary to previous trauma). PECA® Compressive screws can be used for arthrodesis in the midfoot, metatarsal and phalangeal bones of the foot, either alone or in combination with an Airlock® plate.

TRADITIONAL SURGERY




MTP




TALO-NAVICULAR

PERCUTANEOUS SURGERY

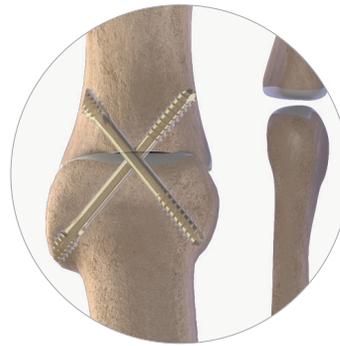

MTP

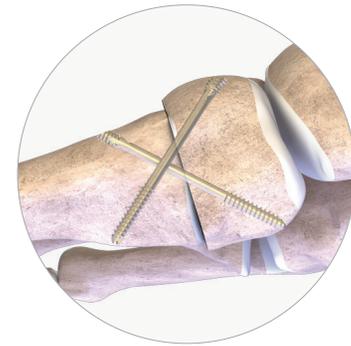

LAPIDUS



MTP ARTHRODESIS POST-OP




MTP ARTHRODESIS POST-OP

KEY FEATURES

Peca® Compressive Implants are intended for the fixation of osteotomy and arthrodesis of the foot. The cannulated, titanium alloy TA6V ELI implant is equipped with a beveled head and partially threaded cutting flutes allow compression by seating flush with the near cortex, preserving soft tissue structures.



BEVELED DESIGN

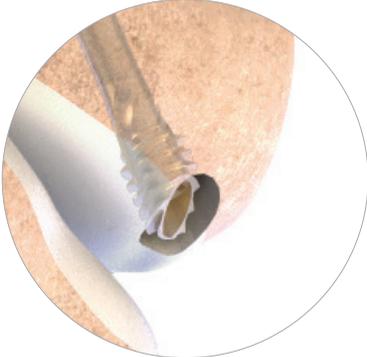

BEVELED HEAD

Cortical anchorage maximization and soft tissues preservation

VS




STANDARD HEAD

Compromise between burial of the head and cortical anchorage

EXACT-2T RECESS

Allows exact beveled implant head positioning. The driver is designed to achieve exact positioning of the screwdriver tip in one position only.

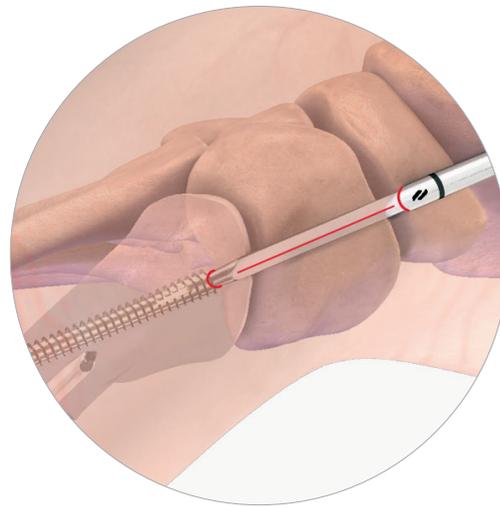

EXACT2-T RECESS



SPECIFIC

Easy indexing of the Exact2-T screwdriver tip.

UNIVERSAL


Possible removal with standard instrumentation.




DEDICATED INSTRUMENTATION

Ø 4mm PECA® Compressive implant can be used:

- With PECA® set, with the optional addition of a countersink;
- With AIRLOCK® set, with the addition of PECA® Compressive instruments.

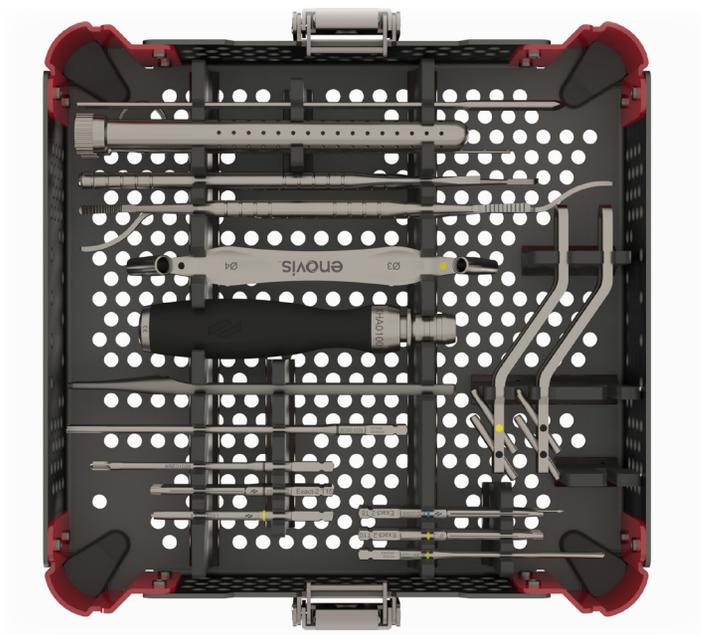

PECA® SET

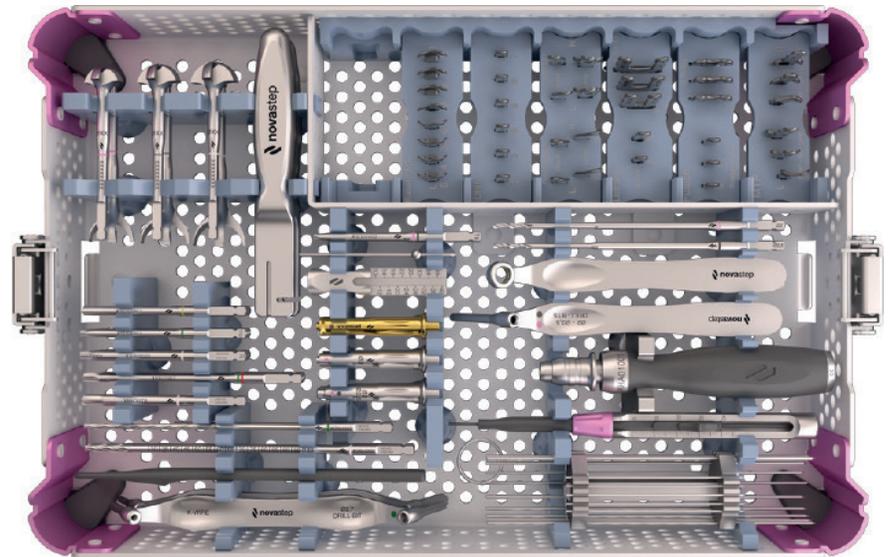

AIRLOCK® SET

REFERENCES

PECA® COMPRESSIVE IMPLANTS

LENGTH	REFERENCE
18mm	PS060118
20mm	PS060120
22mm	PS060122
24mm	PS060124
26mm	PS060126
28mm	PS060128
30mm	PS060130
32mm	PS060132
34mm	PS060134
36mm	PS060136
38mm	PS060138
40mm	PS060140
42mm	PS060142
44mm	PS060144
46mm	PS060146
48mm	PS060148
50mm	PS060150
55mm	PS060155
60mm	PS060160

PECA® COMPRESSIVE SPECIFIC INSTRUMENTS

REFERENCE	DESCRIPTION
CKW01003 ⁽¹⁾	K-WIRE Ø1.6 LG 150 TR/RD
XDB01030	AO DRILL BIT Ø3.2
XSD09002	EXACT2-T15 AO SCREWDRIVER TIP
XGA01009	RULER LG 150
XRE01026 ⁽²⁾	COUNTERSINK Ø3.7

⁽¹⁾ K-wire supplied separately

⁽²⁾ Optional

enovis.

T +33 (0) 2 99 33 86 50 F +33 (0) 9 70 29 18 95

Legal manufacturer: Novastep® S.A.S
2 Allée Jacques Frimot | 35000 Rennes | France
contact-nsp@enovis.com
www.int.novastep.life

Copyright © 2024 Enovis Foot and Ankle

Carefully read the enclosed Instructions For Use (IFU) and all packaging label information. Devices: Implants: Class IIb-CE1639 / Instruments: Class I / Class Ir-CE1639 / Class IIa-CE1639.

REFERENCE: PECA-C-2.0-L-ED4-01-26-EN